



enterMed

MICRONOMIC 50®

For ENT consulting rooms

The budget-friendly MICRONOMIC 50® is a mobile treatment unit designed to accommodate the maximum number of integrated functions according to specific requirements. Providing maximum storage space, the MICRONOMIC 50® incorporates an articulated arm that can be positioned for convenient accessibility to the handpieces. Each handpiece holder is equipped with an infrared sensor which automatically triggers the integrated functions. The two instrument storage drawers are fitted with dividers and trays that slide from side to side in order to optimise instrument storage capacity. Also included are containers for waste and used instruments with letter-box openings. The MICRONOMIC 50® comes with a durable polyurethane worktop and metallic silver lacquer finished front panels. This unit is particularly suitable for ENT consulting rooms with limited space.



**Compact
ergonomic
design**

**Durable
steel
construction**

**Mobile on
castors**

**Ample
storage
space**

**Automatic
on/off of
functions**



enterMed




Standard integrated functions:

- Fibre optic light source with light cable
- Suction system 35 l/min
- 2 mains power sockets for separate devices

Optional integrated functions:

- Compressed air system c/w 2 liquid 1 powder spray bottles
- Ear syringe with electronic temperature control
- Mirror preheater
- Hot air mirror heater
- Heated telescope holders
- Non-heated telescope holders
- Low voltage power source with switchable hook for Frenzel glasses and examination headlight
- Diagnostic microscope with 3-step magnification changer
- Full HD camera system with integral computer and USB connector in console
- 21" touch-screen HD monitor
- 100 Watt HF generator (mono and bipolar)
- X-ray viewer on a post

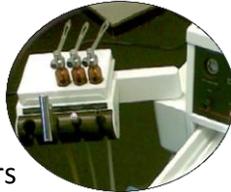

Trays and racks:

A number of trays and racks are included to maximise storage efficiency of both the MICRONOMIC 50® and MICRONOMIC PLUS®.

- Ear speculum rack (extra charge)
- Aluminium instrument trays
- Dividers for drawers

Approx dimensions of MICRONOMIC 50® and MICRONOMIC PLUS®:

Length: 60 cm
 Width: 60 cm
 Height: 90 cm
 Weight: 57 kg



Chassis and side housing in white. Lacquer finished front panels in metallic silver (RAL 9006). Other RAL colours optional. Tough polyurethane top in white (RAL 9016).

MICRONOMIC PLUS®

For storage of instruments

The mobile MICRONOMIC PLUS® is a storage unit particularly suitable for use as an extension to the MICRONOMIC 50®.

The MICRONOMIC PLUS® incorporates:

- Integrated pull-out instrument tray
- 5 large storage drawers with dividers and trays
- Recess in polyurethane worktop for aluminium tray
- Unit lighting



MICRONOMIC PLUS®



enterMed

MICRONOMIC COMBI®

The MICRONOMIC 50® and MICRONOMIC PLUS® can be combined with connecting shelves and a long Marlan worktop to form the MICRONOMIC COMBI®



MICRONOMIC COMBI® with microscope, HD camera system and 100 W HF generator

Approx dimensions of the MICRONOMIC COMBI®:



Height: 90 cm
Length: 160 cm
Depth: 60 cm
Weight: 120 kg (empty)

Mobile on castors

Product code MICRONOMIC 50®: 980554
Product code MICRONOMIC PLUS®: 980562
Product code MICRONOMIC COMBI®: 980566



entermed

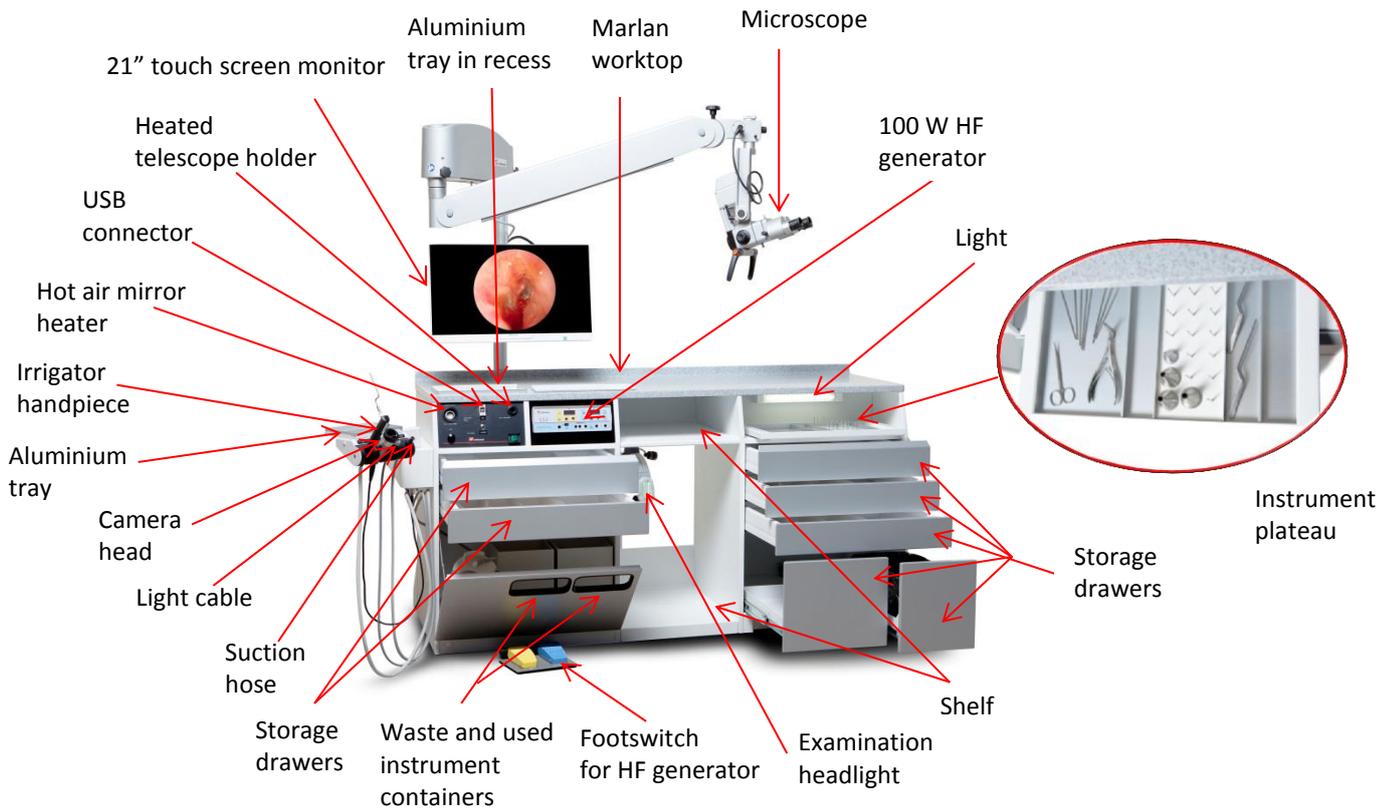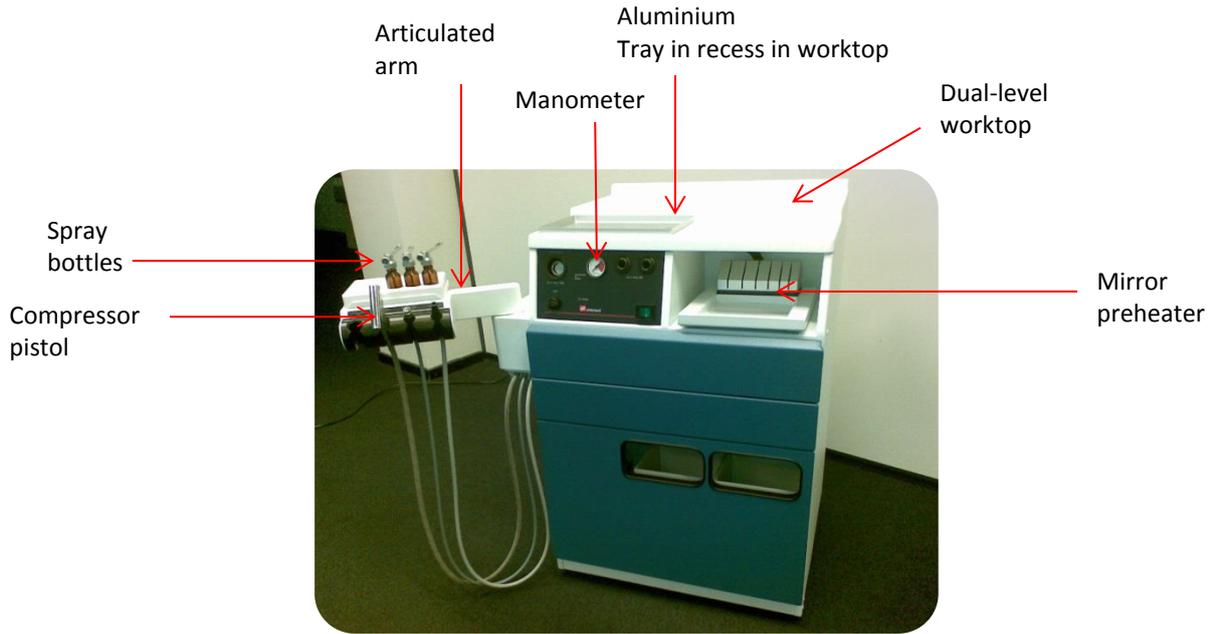

Operates on 230V a/c (±10%) – 50/60 Hz

Entermed B.V.
 Edisonstraat 22
 8912 AW Leeuwarden
 The Netherlands
 Tel: +31 58 298 8888
 Fax: +31 58 298 8768
 Email: info@entermed.com
 Internet: www.entermed.com

We continuously strive to improve and update our products and thus reserve the right to make alterations without prior notice.